



Technology in Practice™



Q200/HE™

Digital Recorder

COMBINATION HOLTER & EVENT RECORDER

ACCURATE AND RELIABLE

CAPTURE ELUSIVE CARDIAC EVENTS

ERGONOMIC DESIGN

Advancing Holter Technology

The Q200/HE™ is the industry's first combination 14-day Holter monitor plus 30-day Event recorder. The Q200/HE interfaces with Holter LX Analysis software to capture elusive cardiac events, deliver fast and accurate analysis, and provide customizable reports. It is the best solution for your practice and your patient, another good combination!

FEATURES & BENEFITS

Combination Recorder: Deployable as a 14-day Holter recorder or a 30-day Event recorder.

Ergonomic Design: Large, easy view LCD screen, soft-keys and a very simple user interface make the unit exceptionally easy to use.

Multiple Language Capability: English, Danish, French, German, Italian, Polish, Portuguese and Spanish.

Low Cost, Shielded Cables: The QRS patient cables provide better signal integrity and noise immunity while the low cost drives down the total cost of ownership.

Accurate and Reliable: Patient data is encoded on a memory card so no information confusion will occur. Reliable digital technology is used throughout recording to enhance service life and lower costs.

Maximize Reimbursements: The Q200/HE allows for "capturing" reimbursement of current procedural Terminology allowed for Holter and Event. At the same time, removable digital storage allows for rapid re-deployment of the unit to another patient.

HOLTER LX ANALYSIS AND LX EVENT

- The Q200/HE interfaces with multiple levels of Holter LX Analysis software, Basic, Enhanced, Enhanced Plus and Pro, each level tailored to address varying needs of our customers. Consult our Holter LX Analysis Software brochure to choose the level that is right for you.
- The Q200/HE interfaces with LX Event software which offers the most comprehensive and flexible way to manage ECG data from looping event recorders.
- Holter LX Analysis and LX Event softwares offer customizable reports, accurate analysis and is available in multiple languages.



WHAT'S INCLUDED

- Q200/HE Recorder
- 5-Lead Patient Cables
- Q200 Operator's Manual
- Hook-Up Guide
- Monitor Pouch & Belt Clip
- Removable Media Memory Card
- Holter Electrodes
- Patient Prep Kit

SPECIFICATIONS

Recording duration: 14 days with alkaline battery

Recording bandwidth: 0.05 – 70 Hz

Resolution: 12-bit, 12.50 μ V/LSB

Sampling rate: 180 samples/second


Options: Language versions available

Dimensions: 8.6 X 6.0 X 2.0 cm

Weight: 99.3 grams, including battery
(70.9 grams without)

Power: 1 AA battery

Storage Capacity: 32 MB – 1 GB SD Flashcard

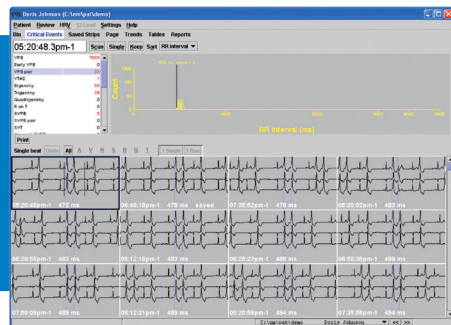

Technology in Practice™

6901 E. Fish Lake Road, Suite 188 | Maple Grove, MN USA 55369
+1.800.465.8408 | +1.763.559.8492 | Fax +1.763.559.2961
info@QRSdiagnostic.com | www.QRSdiagnostic.com



REF 700005-00